

Lower body x-ray shield

312/E-006 (for use on 25 mm or 29 mm high table side rails)

Table mounted shield – with 57 cm detachable top
(compatible with additional optional 57 cm & 27 cm detachable top shields)

This lower body shield substantially reduces exposure to scatter radiation emanating from the table or patient's body. The shield comprises a lightweight and durable structure, with a central friction maintaining pivot link and pivotal arm supporting a flexible lead curtain. A 57 cm wide semi-flexible top shield provides additional protection to the upper body.



Features

- > The shield easily mounts to either side of the table allowing the pivotal arm to always face the head end.
- > A 28 cm wide attachment bracket slides onto a 25 mm high or 29 mm high table side accessory rail (max. 11 mm thick).
- > The central pivot has an automatic wear adjustment facility ensuring accurate joint friction which prevents unwanted curtain movement even when the table is tilted in any direction.

Flexible x-ray protective curtain


- > The curtain forms a large (78 cm wide x 70 cm high) safe zone of occupancy for the Radiologist (or Cardiologist).
- > The curtain forms around an under table x-ray tube without impeding movement.
- > The curtains inner 0.5 mm lead equivalent layer is bonded to its outer covering. As the protective layer is not suspended separately in a cover (bag), it cannot tear or fall away inside unnoticed.
- > The pivotal arm and central pivot can swivel outwards from the table side to form a 50 cm wing.




Optional top shields

- > A 27 cm wide top (312/DS/3.30) can be located onto the attachment bracket.
- > The detachable tops overlap to provide protection to the upper body.
- > These semi-flexible tops stand above the table and can be moved towards the patient as required.
- > Tops can be quickly removed to aid patient access.
- > The x-ray attenuating material used in all tops has a lead equivalent of 0.5 mm.

Weights, dimensions & ordering information

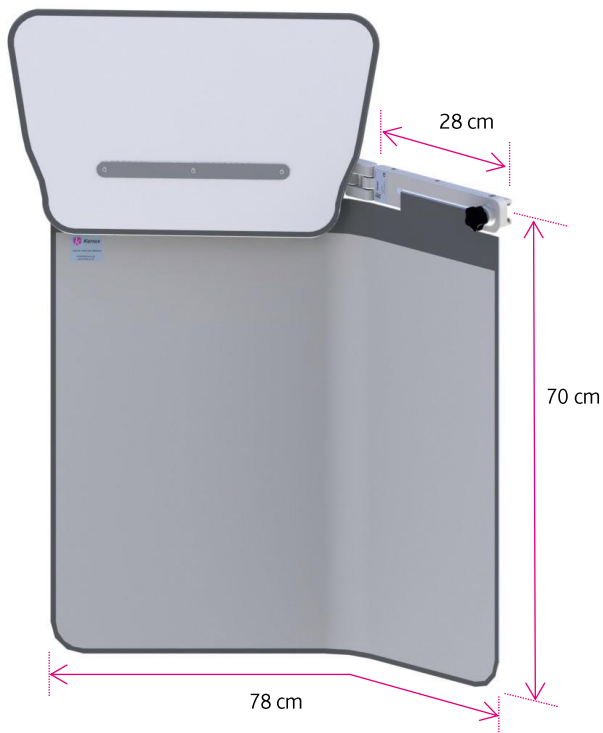

312/E-006
Weight: 6.3 kg

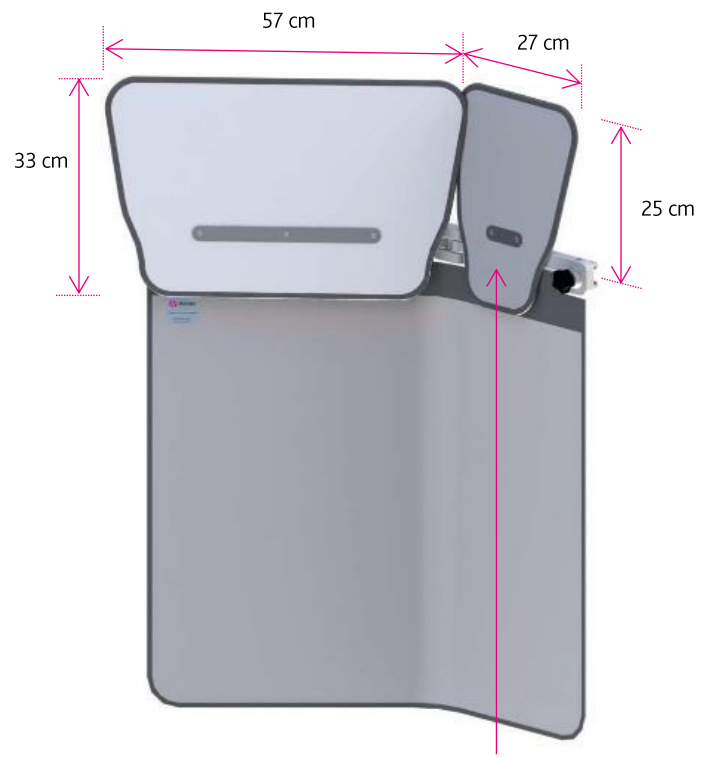

Optional top shield
312/DS/3.30
Weight: 0.8 kg